

Field Service Procedure

GCX Suction Canister Mounting Assembly Anesthesia



**Revision -
12/3/03
SP00299**


Service Procedure

GCX Suction Canister Mtg. Assy

GCX Suction Canister
Mounting Assembly

- 1
- Purpose
- The purpose of this document is to provide instruction for the mounting of the Suction Canister Mounting Assembly Kit onto a 10mm rail (not provided).
- 1.1
- List of Parts and Mounting Hardware
- The following parts are included within this installation kit.

ITEM	QUANTITY	DESCRIPTION
1	1	Canister Mounting Bracket
2	1	Canister Hoop

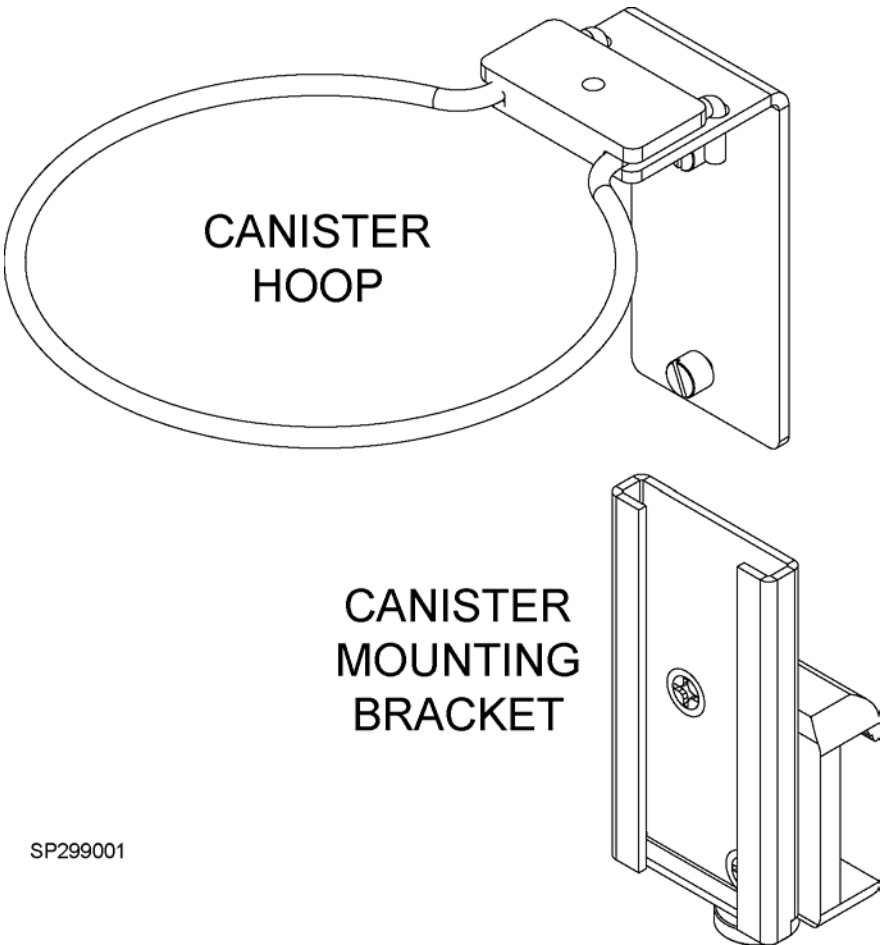

SP299001

Fig. 1 GCX Suction Canister Mounting Assembly


GCX Suction Canister Mtg. Assy

Service Procedure

2 Installation Procedure

2.1 Mounting the Suction Canister Bracket

1. Hook the Canister Mounting Bracket over the 10mm rail (not provided). Tighten the clamping knob (bottom of Mounting Bracket) to secure the Bracket to the rail.

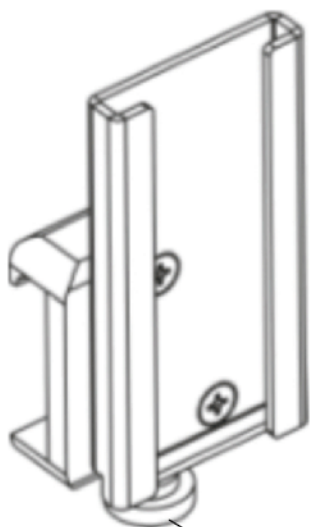

CLAMPING KNOB

SP299002

Fig. 2 Mounting and Suction Canister Bracket

2. Mount the Suction Canister Holder (with mounting detail) in the channel.

Service Procedure

GCX Suction Canister Mtg. Assy

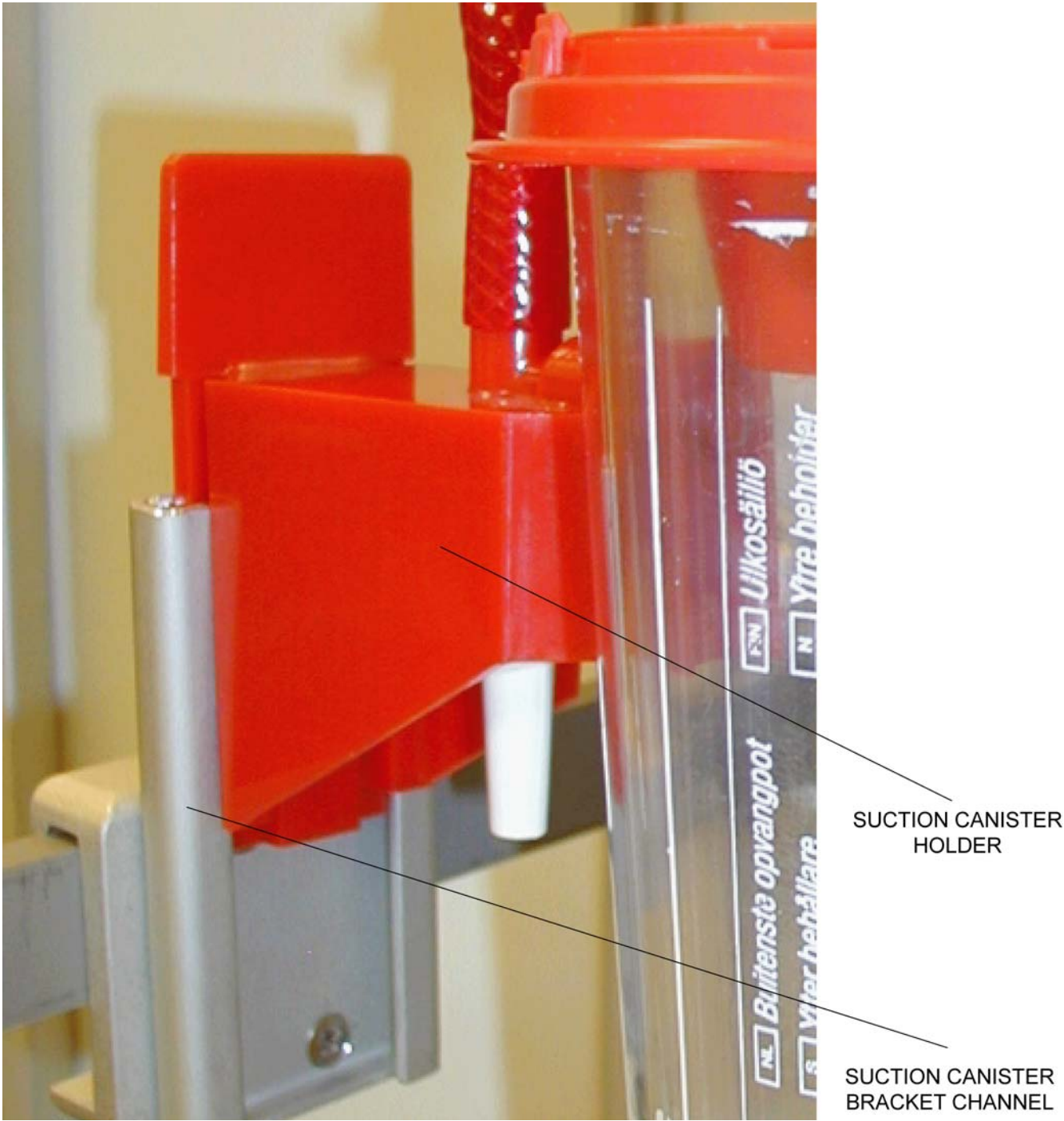


Fig. 3 Mounting the Suction Canister to the Channel

GCX Suction Canister Mtg. Assy

Service Procedure



NOTE:

If the Suction Canister does not provide a mounting detail, refer to Section [2.2 Mounting the Canister Hoop](#).

2.2 Mounting the Canister Hoop

1. Slide the Canister Hoop into the channel and tighten the two nylon screws.
2. Insert the Suction Canister in the Canister Hoop.

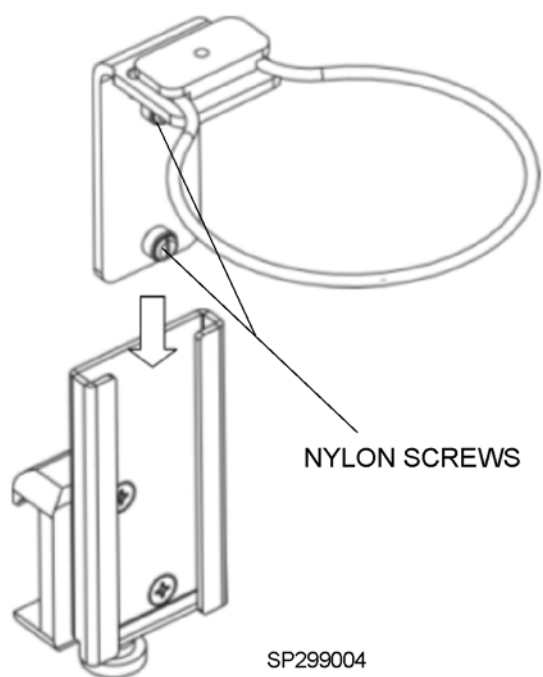


Fig. 4 Mounting the Canister Hoop

In the United States:
Draeger Medical, Inc.
3136 Quarry Road
Telford, PA 18969
USA

Tel: (215) 721-5402
(800) 543-5047
Fax: (215) 721-5785

GCX Suction Canister Mounting Assembly en
SP00299

Subject to change without notice.
Will not be replaced in the event of modifications.

In all other countries:
Dräger Medical AG & Co. KGaA
Moislinger Allee 53 c 55
D-23542 Lütbeck
Germany

Tel: (++49) (0) 1805-3723437
Fax: (++49) 451/882 - 3779